



Cleaning and maintenance of your ThinTube/Tip

Please Note: Do not remove LifeTip from the tube


ThinTubes will need to be replaced every 4-5 months. Spare tubing will issued from the Audiology Department to save patients from returning for each replacement. For removal of any blockages, remove the tubing from the hearing aid and insert the blue cleaning wire in to the end of the ThinTube which attaches to the hearing aid and fed through to the LifeTip. Insertion into the LifeTip whilst attached to the hearing aid could push a blockage further into the tubing and potentially damage the hearing aid.

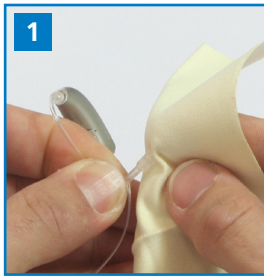

Carefully wipe the LifeTip with a soft damp cloth.

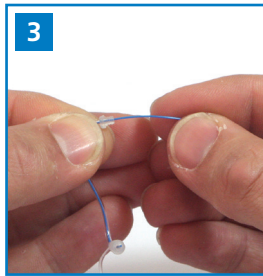


Remove the LifeTip from the ThinTube. Insert the cleaning wire into the ThinTube from the end that has been unscrewed from the hearing aid.

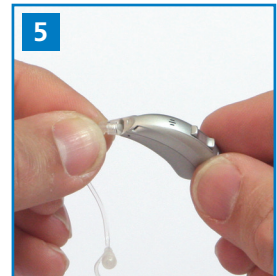

Secure the LifeTip to the ThinTube. Carefully screw the ThinTube to the hearing aid to reconnect.



Carefully unscrew the ThinTube from the hearing aid.

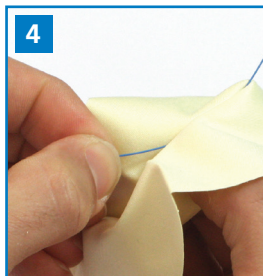

Wipe the cleaning wire with a soft cloth.