

Ear mould hearing aids

Putting your hearing aid in

- 1** Raise the aid to your ear as shown. Tuck the top point of the mould into the fold of your ear.
- 2** Put the bottom half of the mould down into your ear.
- 3** Press the mould piece firmly into your ear.
- 4** You may need to use your other hand to stretch your ear to ease fitting.
- 5** Hook the aid over the top of your ear. Turn the hearing aid on.



Inserting Your Ear Mould



1 Hold the ear mould by the back curved part like a cup. Making sure the hooked part is in an upwards direction.

2 Insert this part under the fold at the top of the ear and at the same time insert the round bulbous portion with holes in down your ear canal.



3 Loosen your grip of the curved portion and fold it backward and press firmly into the bowl of your ear.



Cleaning and maintenance of your ear mould

- 1** Carefully remove the ear mould from the hearing aid.
- 2** Carefully wash the ear mould under a running tap.
- 3** Carefully dry the ear mould either with a soft cloth or a gentle shake. To dry the tubing either use an air blower or leave to dry naturally. Please do not reconnect the ear mould until completely dry. Please ask your Audiologist for more information about a cleaning kit which you may purchase from the Audiology Department.
- 4** Carefully reconnect the ear mould to the hearing aid.
- 5** Ensure the ear mould and hearing aid are connected correctly. Align correctly by placing the ear mould curve to the inside curve of the hearing aid.

***Please Note:** When the tubing becomes damaged contact the Audiology Department and they will replace it for you.*

Do not wash your hearing aid; a wipe over with a baby wipe is sufficient. Your ear mould should be serviced every 6-9 months.


Cleaning Your Earmould

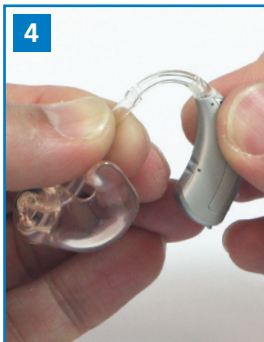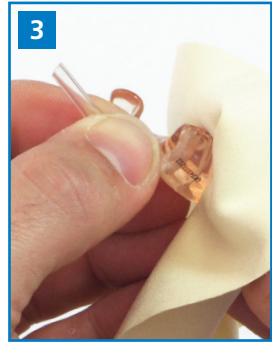


- 1** Remove your ear mould from the hearing aid.
- 2** Wash ear mould in warm soapy water and rinse.
- 3** Dry thoroughly, do not blow into the tube.
- 4** Re-attach the tube to the appropriate hearing aid.

Please Note: If you are unable to clean your hearing aid, please contact your hearing care professional for a service. Use the open access drop-in service for this.



How to retube your ear mould

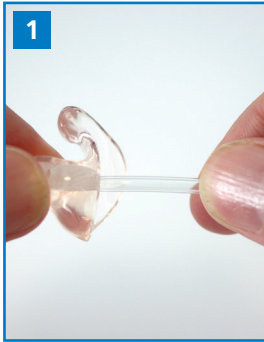

1
Pull old tubing out and keep for later.

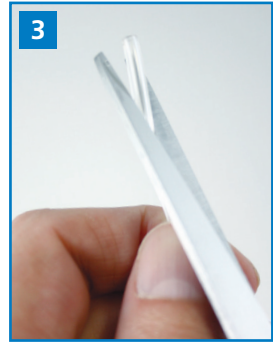


3
Cut end off new tubing into a taper.

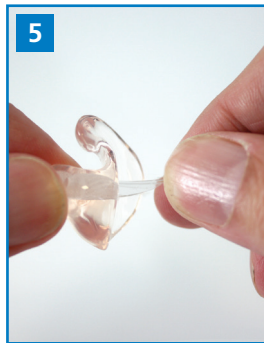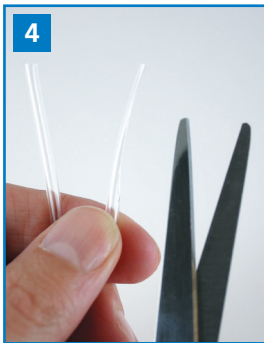


4
5
Push tapered end into ear mould and feed through the hole. If there are two holes in your ear mould, use the top hole.




6
Feed through until bend reaches ear mould and tube is in an upright position.

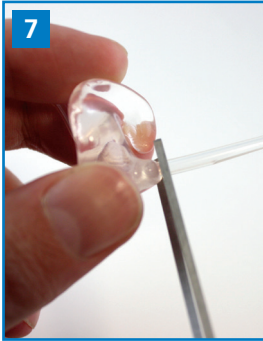

7
Cut tapered end close to the ear mould.

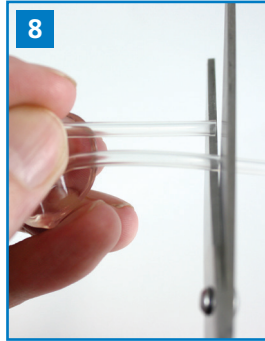


8
Line up with old tubing and cut...

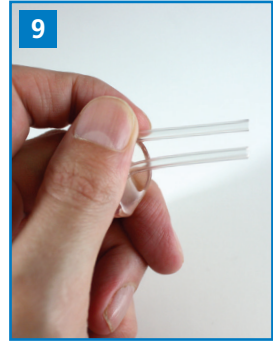

9
...to same length



10
Reattach hearing aid to ear mould.

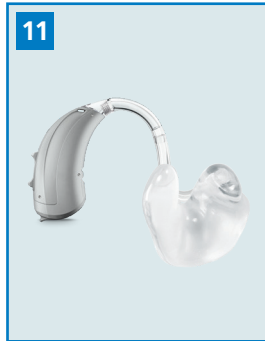

11
Curve of the ear mould should line up with curve of the hearing aid.